




Discriminate Spinal Immobilization *XCollar Field Trial*

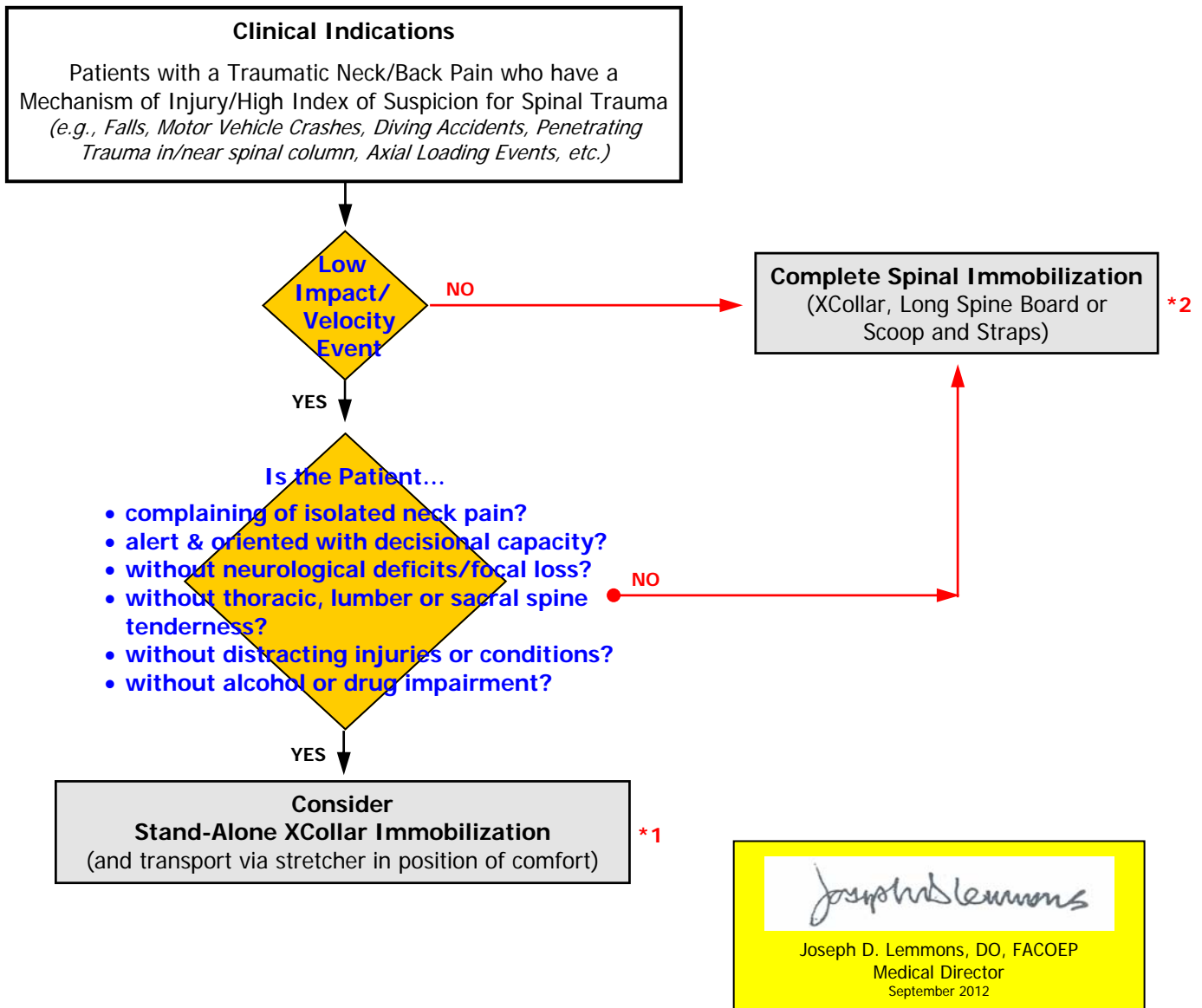

Pearls

*1 XCollar is a unique cervical splinting system designed to splint the Cervical Spine by securing the head to the torso of the patient above C-1 and below C-7; on two points anterior and two points posterior. This splinting system utilizes bi-lateral and vertical adjustments to obtain a customized fit for both the desired circumference and height of large adult and small pediatric patients; working to increase patient safety by minimizing potential Cervical Spinal distraction and injury.

*2 XCollar has an Integrated Head Restraint System (HRS), which replaces the use of head blocks, towels, or tape. KED and HID are not required and can be deferred when using the XCollar.